

# Pharmacy Plus Bluetooth Enabled Refrigerator

## 310L Pharmacy Plus Upright Fridge with Solid Door with Digital Lock

GDMA code PPSR310BT-DLK

SKU code 444411391




### Key features

- ✓ Local monitoring and alarm history via Bluetooth
- ✓ Solid and glass door options
- ✓ Power failure alarm with backup battery powered recording
- ✓ Intelligent fan management system for enhanced performance



#### Bluetooth Connected Controller

Local temperature monitoring and alarm history via Bluetooth connectivity. Power failure alarm with battery back-up, providing continuous temperature monitoring in the event of power outage.

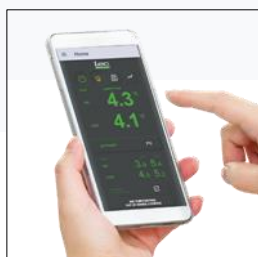


#### Lec Medical Smartphone App

Contactless local monitoring of air and load temperatures via smart phone app. Downloadable temperature report and alarm history.



#### Intelligent Fan Management

Smart fan optimisation aids in temperature stability following door openings, minimising fan activity when door is open to reduce the rate of temperature increase.

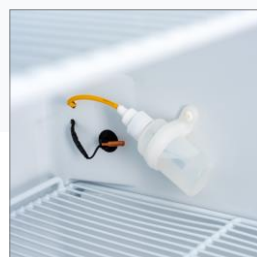

#### Dual Probe Temperature Monitoring

Our dual probe system offers greater temperature accuracy within your appliance, simultaneously measuring both air and product temperatures.



#### Flexible Internal Storage

Multiple shelf locations provides flexible storage solutions.

Specification: 310L Pharmacy Plus Upright Fridge with Solid Door and Digital Lock

Gross Capacity (litres)	310L
Net Capacity (litres)	305L
Unit External Width (mm)	595mm
Unit External Depth (mm)	660mm
Unit External Height (mm)	1560mm
Product Weight (kg)	50.0Kg
Packed Width (mm)	640mm
Packed Depth (mm)	660mm
Packed Height (mm)	1560mm
Packed Weight (kg)	57Kg
Country of Origin	UK
Internal Dimensions (mm) (WxDxH)	515mm x 480mm x 1415mm
Temperature Range	+2°C to +8°C
Colour	White/green
External Finish	Painted Steel
Internal Finish	HIPS*
Type of plastic	SRL800 (HIPS/ABS mix)
Door Finish	Painted Steel / Solid
Door Type	Solid
Door Style	Hinged, Reversible
Door Liner	Flat
Number of Door Trays	n/a
Number of Door Tray Positions	n/a
Door Open Depth (mm)	1150mm
Replacement Door Seal	Yes
Open Door Sensor	Yes
Lock	Yes, 2 Keys
Adjustable Feet (number & position)	2 adjustable (front)
Porthole	Yes
Porthole Size	16 mm
Porthole Location	Rear Bottom
Temperature Display on Panel	Yes
Temperature Alarm	Yes
Remote Alarm Contacts	Yes
Audible Alarms	Yes
Visual Alarms	Yes
Number of Shelves/Drawers	5
Number of Shelf Positions	9
Shelf Size (mm) (WxDxH)	496mm x 405mm x 6mm
Shelf Material	Polythene Coated Wire
Drawer Size (mm) (WxDxH)	n/a
Drawer Material	n/a
Number of Baskets	1
Basket Material	Polythene Coated Wire
Basket Size (mm) (WxDxH)	485mm x 255mm x 171mm
Interior Light	Yes
Light Type	LED
Light Activation	Door Switch / Controller
Fan	Yes - Tangential
Spark Free	No
Defrost setting	Auto (Off Cycle)
Controller Type	Digital BT, External
Battery Back Up/Hours	No
Bluetooth	Yes
NFC	Yes
WIFI	No
Number of Probes	2
Probe Type	Air & Load (Silicone Oil)
Ambient Temperature Range	10°C to 32°C
Refrigerant Type	R600a
CFC Free	Yes
HFC Free	Yes
HCFC Free	Yes
Ammonia Free	Yes
UK 3 Pin Plug	No
Other Plug Types Available	Yes
Voltage Supply	220-240V
Insulation Thickness (mm)	45mm
Insulation Type	High Density Polyurethane Foam
Blowing Agent	Cyclopentane
Anti Microbial Handle	Yes
Wall Brackets	No
Optional Drawer System	No
Optional Digital Lock	No
Other Door Option	Glass PPGR310BT



# WolfLabs

**Pricing on any accessories shown can be found by keying the part number into the search box on our website.**

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

**[www.wolflabs.co.uk](http://www.wolflabs.co.uk)**

**Tel : 01759 301142**

**Fax : 01759 301143**

**[sales@wolflabs.co.uk](mailto:sales@wolflabs.co.uk)**

Please contact us if this literature doesn't answer all your questions.